AlproJet-Series

Our team for cleaning and disinfection of suction systems:
AlproJet-D with AlproJet-W or AlproJet-DD with AlproJet-W

AlproJet-D

Liquid concentrate for the daily disinfection and cleaning of dental suction systems with and without amalgam separator.

AlproJet-DD

Liquid concentrate for the daily disinfection and cleaning of dental suction systems with and without amalgam separator.

AlproJet-W

Liquid concentrate for the weekly acid contra-cleaning of dental suction systems with and without amalgam separator.



ALPRO®
ALPRO MEDICAL GMBH

Cleaning and disinfection of suction systems with AlproJet Series

Recommendation and dosage



1



After each treatment aspirate 1 – 2 cups of cold water and activate the spittoon bowl rinsing for 30 seconds.

With amalgam separator in unit:

Mon	Tue	Wed	Thu	Fri
W	W	D / DD	D / DD	D / DD

With centralized amalgam separator: alternate weekly

Week 1	Week 2	Week 3	Week 4
W	D / DD	W	D / DD

1 Slightly open the screw cap of the dosage chamber.

2 Press AlproJet bottle for dosing.

3 **AlproJet-D**: 2 x 50 ml
AlproJet-DD: 2 x 20 ml
AlproJet-W: 2 x 50 ml

4 Fill AlproJet-D/-DD/-W into AlproJet MixCup.


5 Fill up with drinking water (room temperature) to a 2 litre ready-to-use solution.

6 Attach suction tubes to the adapters of AlproJet MixCup (spray mist aspiration tube to the adapter with the large opening and the saliva aspirator to the adapter with the small opening) and aspirate 1 litre of the ready-to-use solution.

7 Pour in sips the remaining solution (1 litre) into the spittoon bowl and let it take effect.

8 Contact time: 60 minutes at least, preferably overnight/during the weekend.

Please notice:
Rinse suction tubes and spittoon bowl thoroughly with water just **before starting work.**

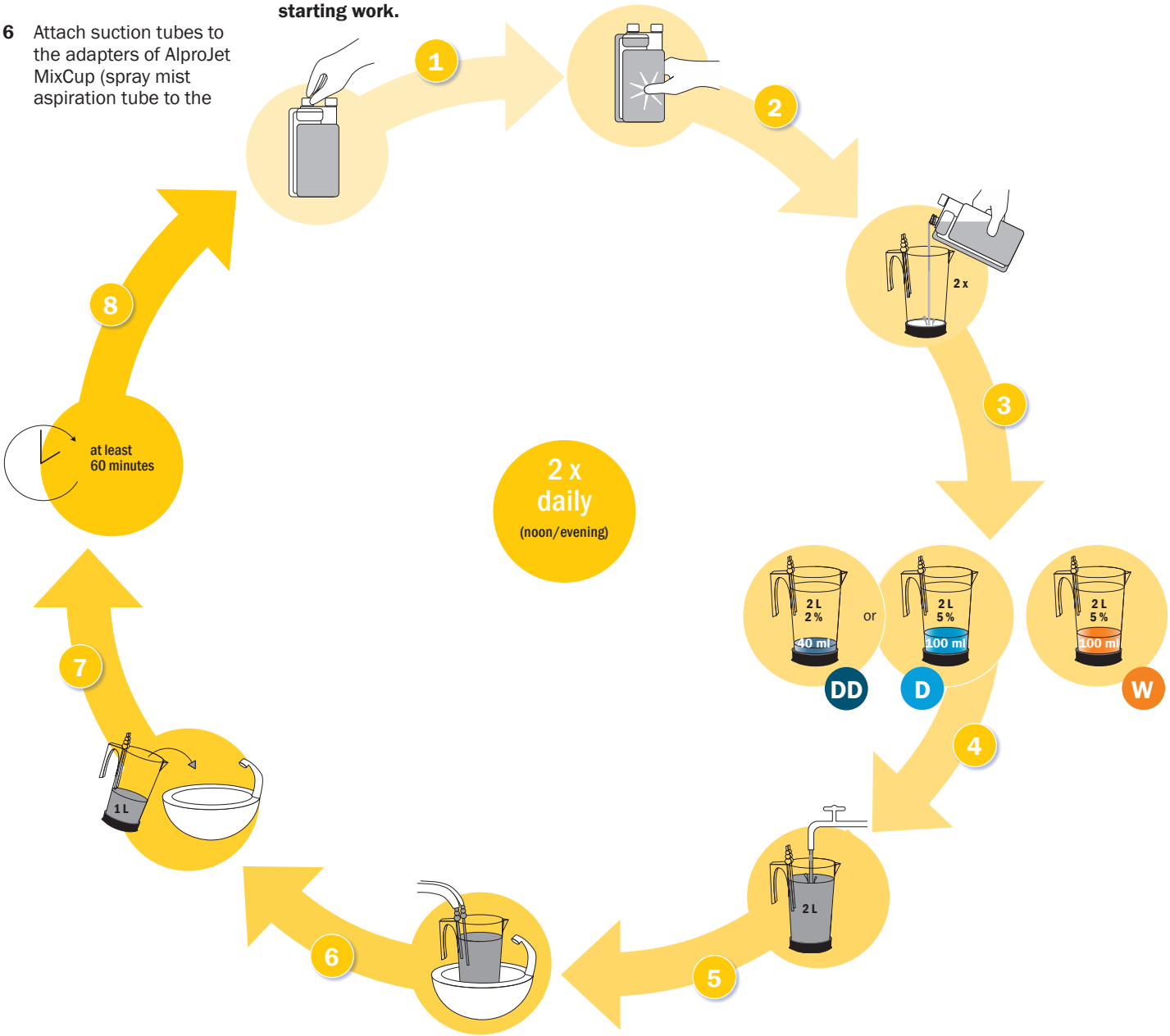

Application of AlproJet-W for surgical work

to prevent deposits, neutralize hydrogen peroxide, reduce foam forming and dissolve coagulated blood



2



After each treatment aspirate 1 – 2 cups of cold water and activate the spittoon bowl rinsing for 30 seconds.

Before surgical work:

1 Slightly open the screw cap of the dosage chamber.

2 Press AlproJet-W bottle for dosing.

3 **AlproJet-W**: 2 x 50 ml

4 Fill AlproJet-W into AlproJet MixCup.

5 Fill up with cold drinking water to a 1 litre ready-to-use solution (10 %).

In case of already coagulated blood in unit:


5a Fill up with drinking water (room temperature) to a 1 litre ready-to-use solution (10 %).

6 Dip suction tubes into AlproJet MixCup and aspirate ready-to-use solution.

Please notice:
Rinse suction tubes and spittoon bowl thoroughly with water just **before starting work.**

Result:

The dissolving process is accelerated, residues become pumpable and thus are removed more easily.

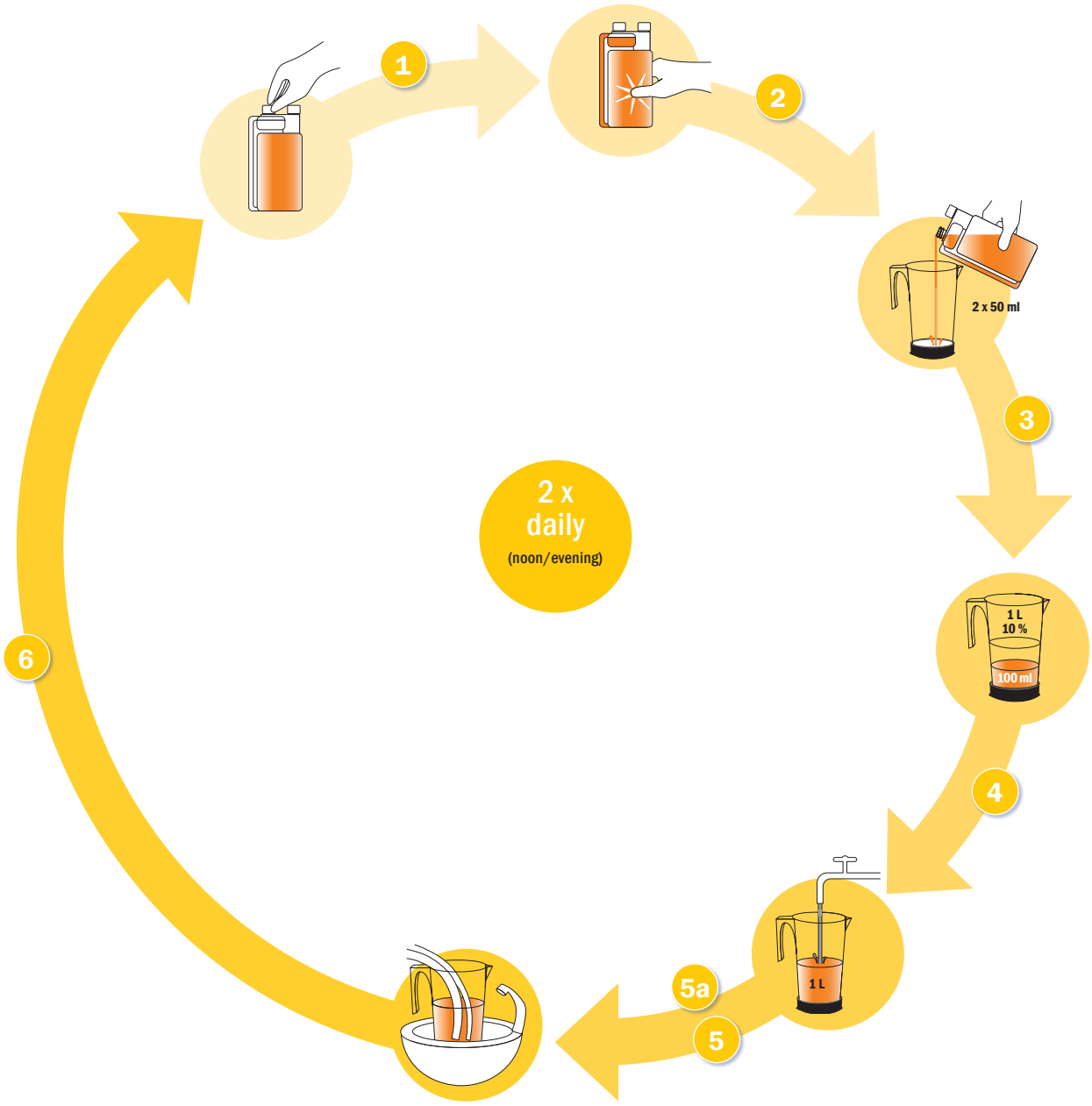

Application of AlproJet-W before prophylaxis treatment

to avoid incrustation, foam forming and suction dropouts, e. g. by powder jet systems



3



After each treatment aspirate 1 – 2 cups of cold water and activate the spittoon bowl rinsing for 30 seconds.

1 Slightly open the screw cap of the dosage chamber.

2 Press AlproJet-W bottle for dosing.

3 **AlproJet-W**: 4 x 50 ml

4 Fill AlproJet-W into AlproJet MixCup.

5 Fill up with cold drinking water to a 2 litre ready-to-use solution (10 %).

6 Attach suction tubes to the adapters of AlproJet MixCup (spray mist aspiration tube to the adapter with the large opening and the saliva aspirator to the adapter with the small opening) and aspirate 1 litre of the ready-to-use solution.

7 Pour in sips the remaining solution (1 litre) into the spittoon bowl.

Please notice:
Rinse suction tubes and spittoon bowl thoroughly with water just **before starting work.**

Result:

During treatment foaming residues are neutralized/eliminated and suction dropouts are prevented.



Basic cleaning with AlproJet-W when changing products

apply 2 x daily for one week



4

1 Slightly open the screw cap of the dosage chamber.

2 Press AlproJet-W bottle for dosing.

3 **AlproJet-W**: 4 x 50 ml

4 Fill AlproJet-W into AlproJet MixCup.

5 Fill up with drinking water (room temperature) to a 2 litre ready-to-use solution (10 %).

6 Attach suction tubes to the adapters of AlproJet MixCup (spray mist suction tube to the adapter with the large opening and the saliva aspirator to the adapter with the small opening) and aspirate 1 litre of the ready-to-use solution.


7 Pour in sips the remaining solution (1 litre) into the spittoon bowl and let it take effect.

8 Repeat point 3 – 7 (in total 2 x).

9 Let it take effect for 1 – 2 hours and rinse with 40 – 50 °C warm water.

Result:

The dissolving process is accelerated, residues become pumpable and thus are removed more easily.

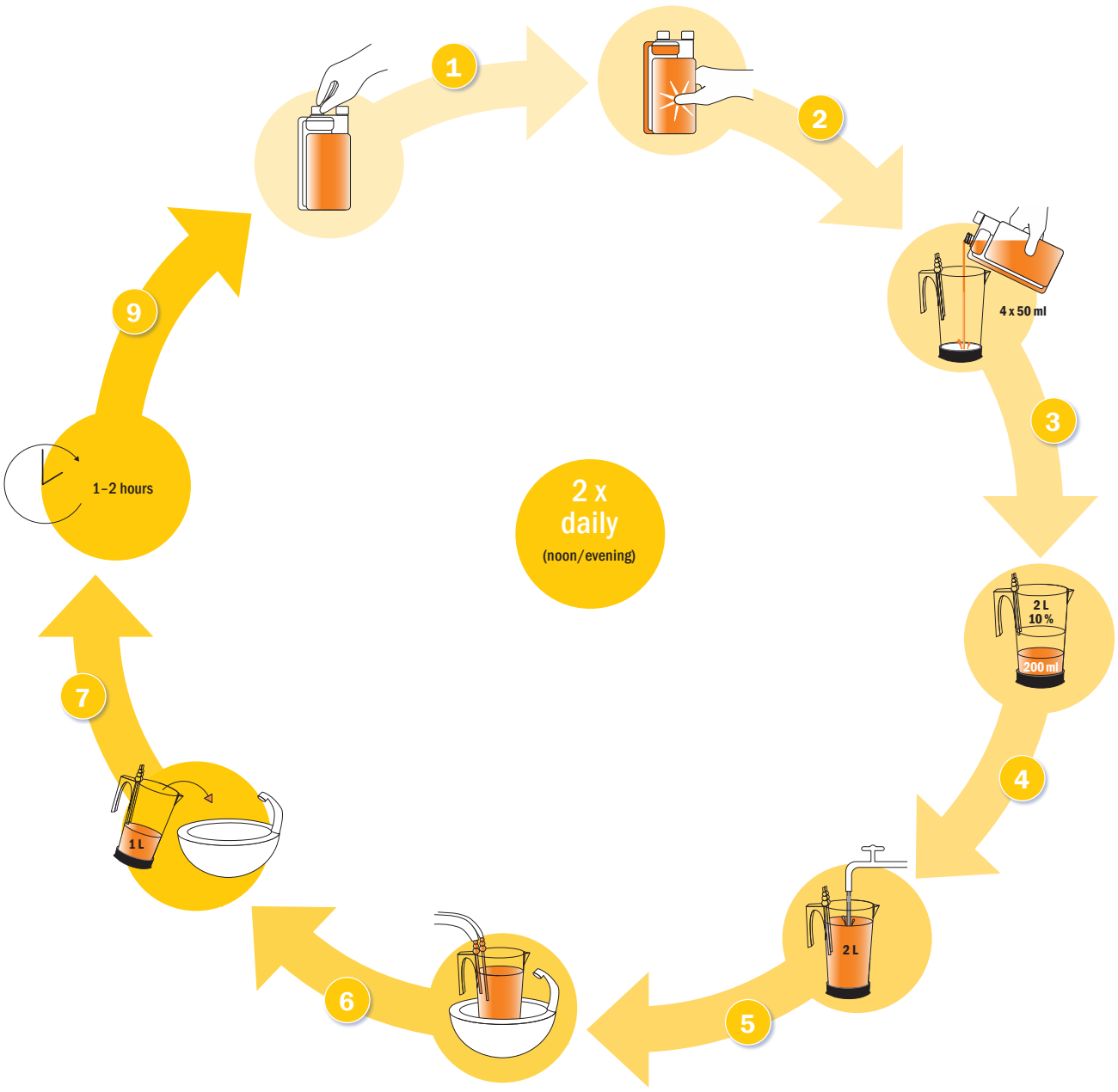



Our solutions for your entire practice hygiene





Our certified and eco-friendly products cover the entire spectrum of infection prevention in clinical practice — from hand, surface and instrument disinfection to the removal and prevention of biofilm in the procedural waterlines of dental treatment units. ALPRO is one of the world's leading manufacturers in the cleaning, disinfection, and maintenance of dental suction systems. Discover our extensive product range and services at:
www.alpro-medical.de

ALPRO MEDICAL GMBH

ALPRO MEDICAL GMBH
Mooswiesenstr. 9
78112 St. Georgen
Deutschland

Tel. +49 7725 9392-0
Fax +49 7725 9392-91
www.alpro-medical.de
info@alpro-medical.de

REF 80350031_05_20260427_GB

Stiftung Zentrale Stelle VERPACKUNGSREGISTER
registered as manufacturer of
system participating obliged packaging