INSTRUCTIONS FOR USE


PlastiSept eco

Surface Cleaning and Disinfection of Medical Devices

Class IIb medical device solution - alcohol-free, ready to use

REF 3063-E · 5 L canister · REF 4645 · 500 mL bottle · ARTG 164389

ALPRO-AU-IFU-PLASTISEPT-ECO · Version 1.1 · Updated October 2026

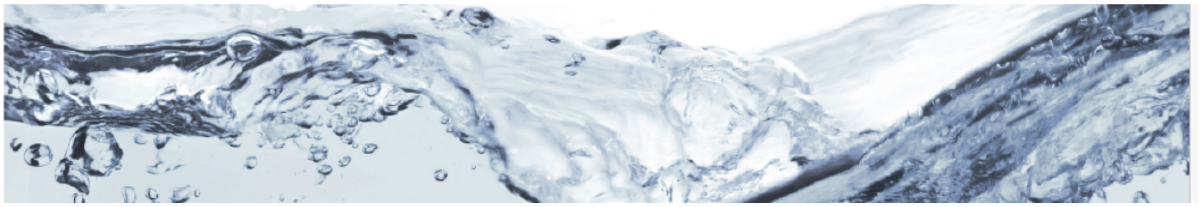

Alcohol-free cleaning and disinfection for sensitive medical-device surfaces

PlastiSept eco is an alcohol-free Class IIb medical device disinfectant suitable for surfaces of dental units, X-ray machines, suspension tables and other non-critical medical devices where disinfection is required and an alcohol-free product is preferred. It is ready to use and formulated for sensitive materials.

Instructions for Use

1. **Ready to use**
   Use PlastiSept eco undiluted. Do not mix with detergents or other chemicals.
2. **Apply to the cleaning cloth - never directly to the device**
   Do not spray or dispense PlastiSept eco directly onto the surface of a medical device. Spray or dispense the product onto a suitable cleaning cloth first, then apply the moistened cloth to the device surface.
3. **Moisten thoroughly without overwetting**
   Wipe the medical-device surface so it is thoroughly moistened, while avoiding excessive wetting or pooling of liquid.
4. **Observe the efficacy contact time**
   Keep the surface visibly moist for the required efficacy period: yeasticidal - 1 minute; bactericidal including MRSA - 5 minutes; virucidal against Vaccinia virus as representative of enveloped viruses - 5 minutes. Allow the surface to remain until the required efficacy time has been reached.
5. **Allow to dry naturally**
   After the required contact time, allow the treated surface to dry. Do not wipe the surface afterwards with a dry cloth.

Use with ALPRO Maxiwipes (250 dry wipes)

When using ALPRO Maxiwipes 250 dry wipes, pour 1.5 L of PlastiSept eco into the wipe container. Allow the wipes to soak for at least 30 minutes before first use so the liquid is distributed evenly throughout the roll.

Weekly surface-care recommendation

Once a week, wipe equipment regularly disinfected with PlastiSept eco using a mild soap solution and water, then wipe with clean water as appropriate for the device. This helps return the surface to a more neutral condition and remove accumulated disinfectant residues.

PlastiSept eco is a water-based disinfectant. During repeated use, the water evaporates while active ingredients can remain on the treated surface. For this reason, weekly residue removal is highly recommended, subject to the medical-device manufacturer's cleaning and compatibility instructions.

Composition and compatibility

Dodecyldipropylenetriamine 0.2-1%, N-C12-14-alkyltrimethylenediamine 0.2-1%, deionised water, fragrances and emulsifying agent. Free of alcohol, aldehydes, QAC and halogen. Always follow the device manufacturer's cleaning and material-compatibility instructions.