INSTRUCTIONS FOR USE


PlastiSept eco Wipes

Surface Cleaning and Disinfection of Medical Devices

Class IIb medical device solution - alcohol-free, ready to use

REF 4625 · 2 tubs + 2 refill bags · REF 4624 · 12 refill bags · ARTG 293206

ALPRO-AU-IFU-PLASTISEPT-ECO-WIPES · Version 1.2 · Updated October 2026

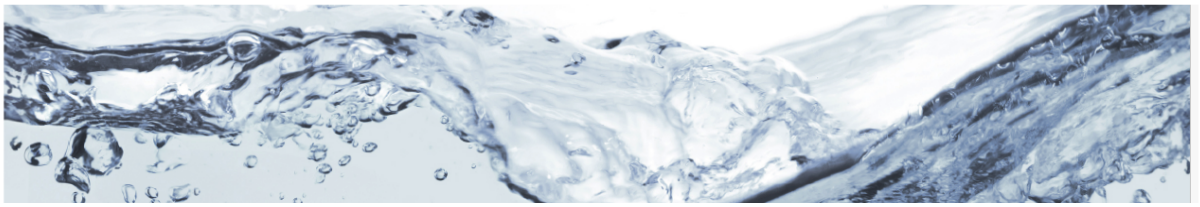

Alcohol-free cleaning and disinfection for sensitive medical-device surfaces

PlastiSept eco Wipes are ready-to-use alcohol-free wipes for pre-cleaning and disinfection of compatible medical-device surfaces. They are suitable for dental units, X-ray equipment, suspension tables and other non-critical medical devices where disinfection is required and an alcohol-free product is preferred.

Instructions for Use

1. **Ready to use**
   Remove one pre-moistened PlastiSept eco Wipe from the dispenser. No dilution or additional solution is required.
2. **Wipe the medical-device surface**
   Apply the wipe directly to the surface and distribute the solution evenly. Thoroughly moisten the surface, but avoid overwetting, pooling or excess liquid around openings and controls.
3. **Observe the efficacy contact time**
   Keep the treated surface visibly moist for the required efficacy period: yeasticidal - 1 minute; bactericidal including MRSA - 5 minutes; virucidal against Vaccinia virus as representative of enveloped viruses - 5 minutes.
4. **Allow to dry naturally**
   After the required contact time has been reached, allow the treated surface to dry naturally. Do not wipe the surface afterwards with a dry cloth.
5. **Close the dispenser after use**
   Close the lid tightly after each wipe is removed to prevent drying. After opening, the wipes remain usable for at least 28 days when the lid is kept tightly closed.

Weekly surface-care recommendation

Once a week, wipe equipment regularly disinfected with PlastiSept eco products using a mild soap solution and water, then wipe with clean water as appropriate for the device. PlastiSept eco is water-based; with repeated disinfection, water evaporates while active ingredients may remain on the surface. Weekly residue removal is therefore recommended to help return the surface to a more neutral condition.

Wipe characteristics

- **Large format**
  21 x 26 cm for greater surface coverage per wipe.
- **Viscose microfibre**
  50 g/m² 'Cotton Touch' material for excellent liquid absorption and even solution distribution.
- **Material compatibility**
  Alcohol-free, aldehyde-free, QAC-free and halogen-free chemistry for sensitive medical-device surfaces.
- **Environmental profile**
  Viscose wipe material offers a more environmentally compatible alternative to many generic polyester wipes.

Soaking solution and compatibility

Dodecyldipropylenetriamine 0.2-1%, N-C12-14 alkyltrimethylenediamines 0.2-1%, deionised water, fragrances and emulsifying agent. Always follow the medical-device manufacturer's cleaning and compatibility instructions.